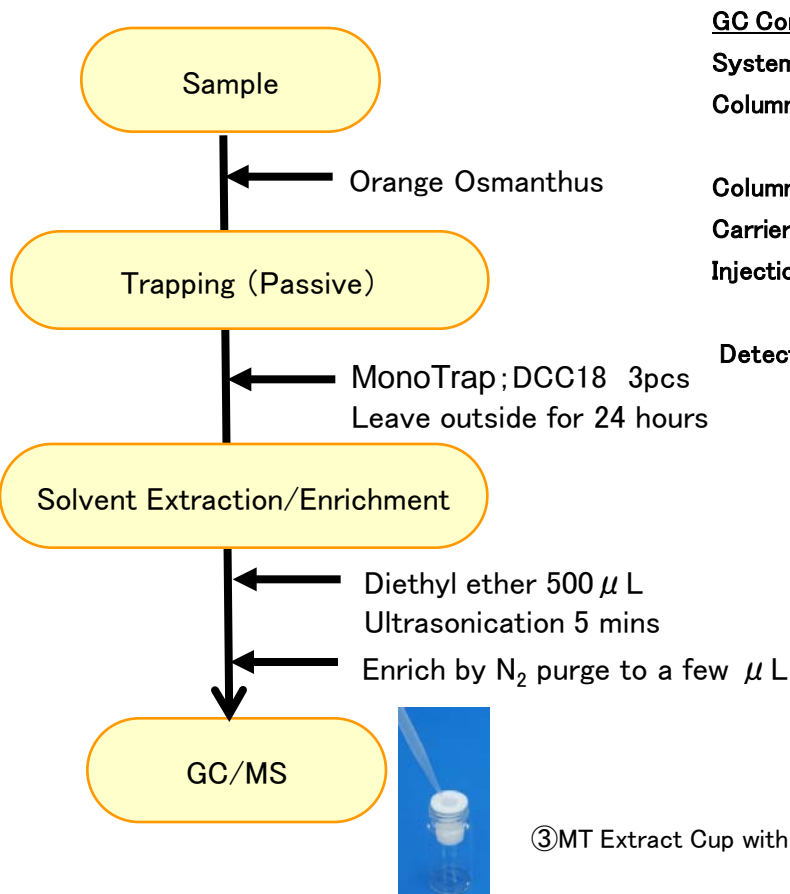


## Protocol

**GC Condition**

**System** : SHIMADZU GC-2010、GCMS-QP2010 Plus

**Column** : InertCap Pure-WAX (Cat.1010-68142)  
0.25mmI.D. × 30m df=0.25  $\mu$  m

**Column Temp** : 40°C(5min)→4°C/min→250°C(5min)

**Carrier Gas** : He 120kPa

**Injection** : Split 1:10, 1  $\mu$  L  
250°C

**Detection** : MS Scan (m/z;35-450)

